





SENOGRAPHE CRYSTAL

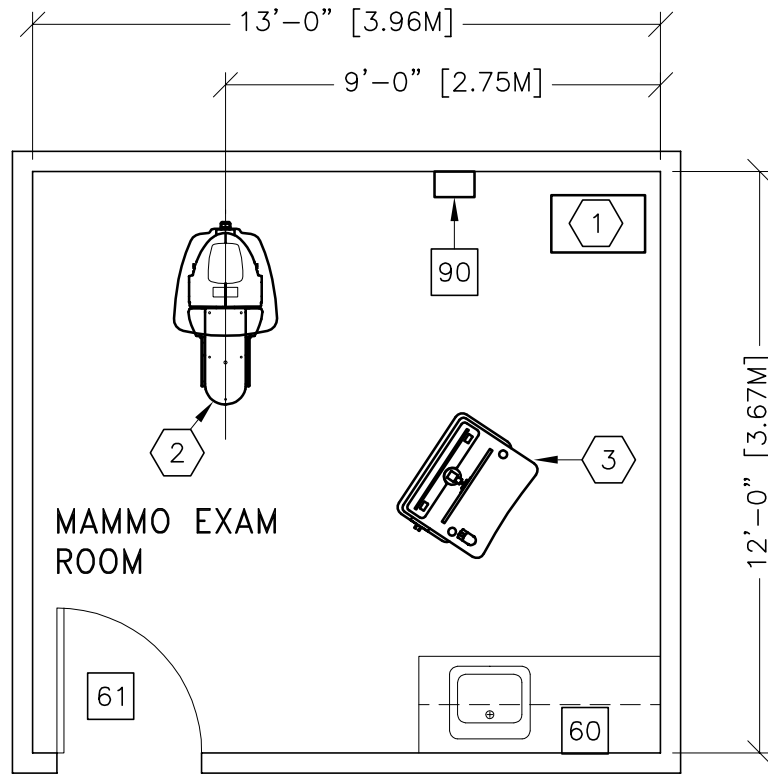

ROOM DIMENSIONS	LENGTH & WIDTH	CEILING HEIGHT
MINIMUM:	9'-11" x 8'-6" [3000MM x 2600MM]	7'-10.5" [2.4M]

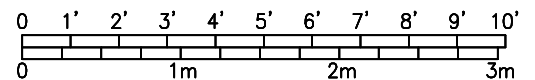

ROOM PLAN LOCATED ON PAGE 9-36

ITEM ⬡	QUANTITY	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
①	1	GENERATOR	82 kg 181 lbs	
②	1	GANTRY (WITHOUT DETECTOR) INCLUDES WATTS OF GENERATOR AND CONTROL STATION	301 kg 664 lbs	441 W 1505 btu
③	1	CONTROL STATION	121 kg 267 lbs	
***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS *****				
60	COUNTER TOP WITH SINK, BASE AND WALL CABINETS			
61	MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 750MM W X 2135MM H CONTINGENT ON A 2500MM CORRIDOR WIDTH			
90	MAIN DISCONNECT CONTROL, WALL MOUNTED GEMS CAT. NO. E4502B, 29kg. SEE DETAIL R4502B.			

ENVIRONMENT

Ambient Operation Temperature: 50-95° F (10-35° C)
 Allowable Temperature Change: N/A
 Humidity: 30-75 percent non condensing
 Allowable Humidity Change: N/A
 Max - Altitude: 9840 ft. (3000M) above sea level