

# GoldSeal Infinia II

Catalog - S9600PT

**The GoldSeal\* refurbished Infinia\* II is a premium, high performance, all purpose dual detector imaging system, with extra large rectangular FOV detectors positioned in flexible geometries.**

## Features Included

Open, path-through flexible geometry gantry

- Slip-ring gantry enables continuous multiple tomographic rotations in the same direction
- Externally mounted detectors, with flexible positioning for all major studies
- Adjustable scanning radius for both 90 degree and 180 degree SPECT scans

Elite digital dual detectors

- Crystal thickness – 3/8 in./9.5 mm
- 59 high quantum efficiency PMT's - each coupled with one ADC per PMT
- Extra large, rectangular FOV: 540 mm x 400 mm (21.25 in. x 15.75 in.)

Dual-axis cantilever imaging table

- Low attenuation carbon fiber tabletop
- Easy patient handling, comfortable tabletop
- Obstruction-free floor installation enables table pivoting-out at a 45 degree angle
- Automatic home positioning for common imaging procedures

Infinia acquisition station

## Configuration

- Infinia II refurbished
- 17" LCD acquisition monitor

## Available options

- Collimators (Preowned LEHR, other Preowned collimators limited availability)
- Table leg extenders
- Pediascan
- Axial head holder
- 10KVA UPS
- IVY R-wave trigger
- Xeleris 4 workstation
- Xeleris 4 software options

## Available Training

- Onsite clinical applications
- TIP Virtual Assist remote training

## Warranty


Includes full 1 year warranty.

## Contact Us

GoldSeal systems are quoted subject to availability. To confirm current availability of systems or additional options in your region, contact your GE Healthcare representative.

Your access to quality medical imaging is our priority. Get the most for your budget by choosing GoldSeal today. If it meets our standards, we know it will meet yours.

Learn more at [www.gehealthcare.com/goldseal](http://www.gehealthcare.com/goldseal)

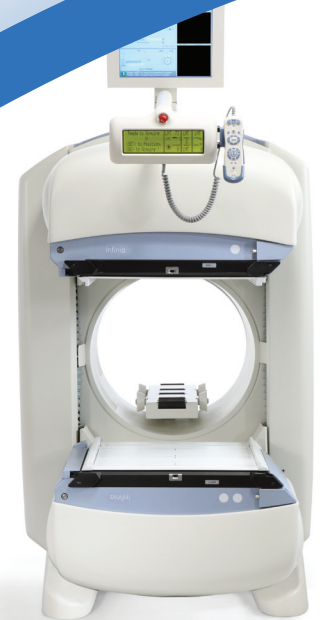

©2016 General Electric Company – All rights reserved.

General Electric Company reserves the right to make changes in specification and features shown herein, or discontinue the product described at any time without notice or obligation. GE and GE Monogram are trademarks of General Electric Company. GE Healthcare, a division of General Electric Company.

\*Trademark of General Electric Company

#### **Good Refurbishment Practice**

The GoldSeal program is consistent with Good Refurbishment Practice (GRP) established by leading industry trade associations such as MITA (Medical Imaging & Technology Alliance) a division of NEMA (National Electrical Manufacturers Association), COCIR (European Coordination Committee of the Radiological, Electromedical and Healthcare IT Industry) and JIRA (Japan Industries Association of Radiological Systems).

## **About GE Healthcare**

GE Healthcare provides transformational medical technologies and services that are shaping a new age of patient care. Our broad expertise in medical imaging and information technologies, medical diagnostics, patient monitoring systems, drug discovery, biopharmaceutical manufacturing technologies, performance improvement and performance solutions services help our customers to deliver better care to more people around the world at a lower cost. In addition, we partner with healthcare leaders, striving to leverage the global policy change necessary to implement a successful shift to sustainable healthcare systems.

Our “healthymagination” vision for the future invites the world to join us on our journey as we continuously develop innovations focused on reducing costs, increasing access and improving quality around the world. Headquartered in the United Kingdom, GE Healthcare is a unit of General Electric Company (NYSE: GE). Worldwide, GE Healthcare employees are committed to serving healthcare professionals and their patients in more than 100 countries. For more information about GE Healthcare, visit our website at [www.gehealthcare.com](http://www.gehealthcare.com)

GE Healthcare  
3114 N. Grandview Blvd.  
Waukesha, WI 53188  
USA  
[www.gehealthcare.com](http://www.gehealthcare.com)



imagination at work