





INFINIA II SYSTEM W/ HAWKEYE OPTION

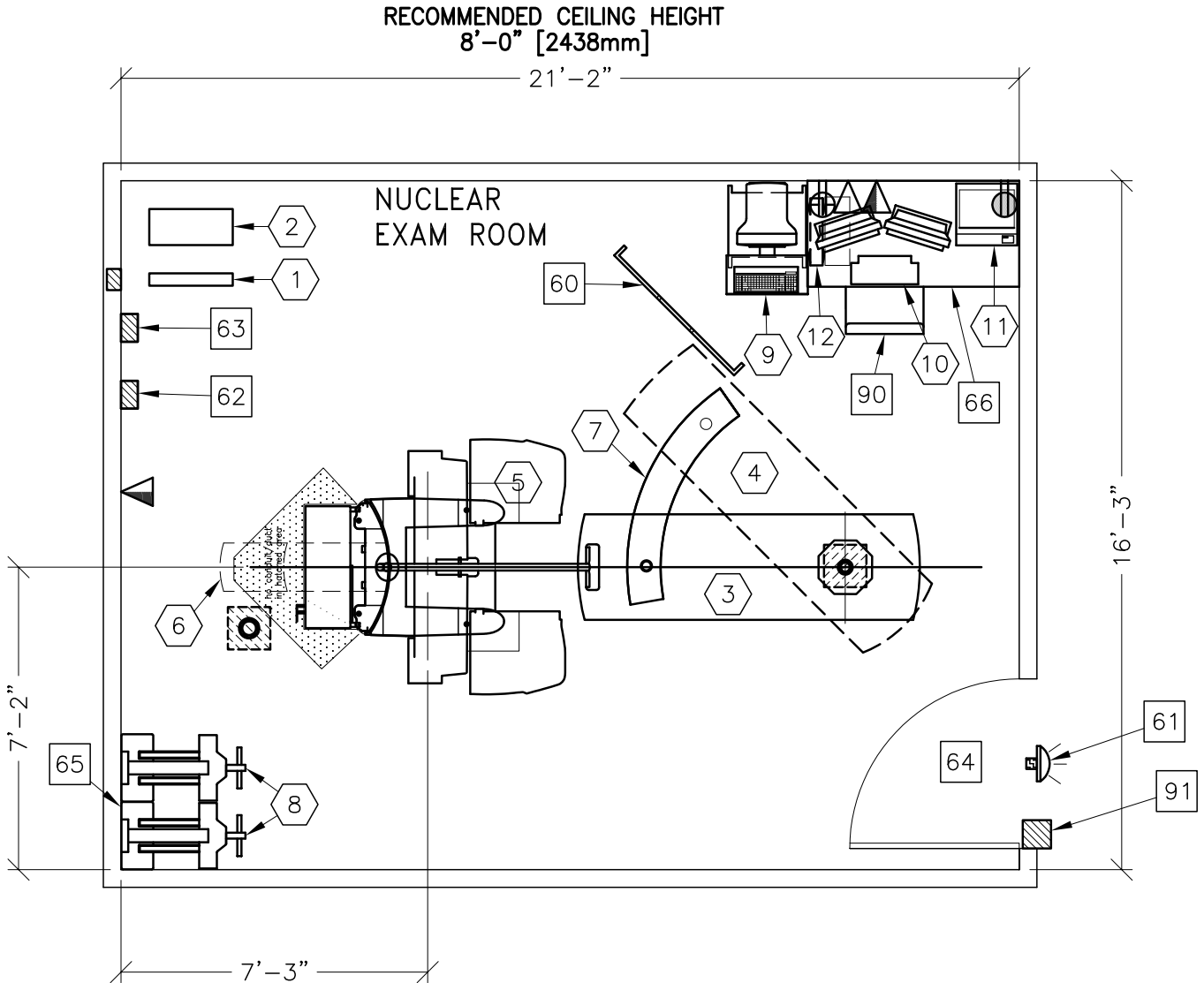

ROOM DIMENSIONS	LENGTH & WIDTH	CEILING HEIGHT
MINIMUM:	Refer to Pre-Installation Manual	

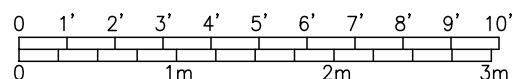

SEE APPENDIX "A" FOR POWER SPECIFICATIONS

ROOM PLAN LOCATED ON PAGE 7-64

ITEM ①	QUANTITY	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
①	2	COLLIMATOR STORAGE CART	1058 lbs	
②	1	INFINIA II IMAGING SYSTEM GANTRY WITH HAWKEYE OPTION	6263 lbs	8600 btu
③	1	IMAGING TABLE	881 lbs	341 btu
④	1	LIMIT OF TABLE TRAVEL		
⑤	1	TABLE SWING FOR COLLIMATOR EXCHANGE		
⑥	1	R-WAVE TRIGGER UNIT	19 lbs	170 btu
⑦	1	XELERIS WORKSTATION	55 lbs	
⑧	1	UPS CABINET (OPTIONAL)	639 lbs	3003 btu
⑨	1	ACQUISITION MOBILE CART INCLUSIVE OF MONITOR AND KEYBOARD	180 lbs	1604 btu
***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS *****				
⑥0		OPTIONAL WALL PROTECTION FROM COLLIMATOR CART		
⑥1		MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 45 IN. W X 80 IN. H [1143MM X 2032MM], CONTINGENT ON A 84 IN. [2134MM] CORRIDOR WIDTH.		
⑥2		OPERATORS CHAIR		
⑥3		CART		
⑥4		TABLE		
⑥5		X-RAY ON WARNING LIGHT AVAILABLE FROM GE SUPPLY CALL: 800-200-9760 EXTENSION 3985 GE CAT. NO. WXIABW-OF-XIU		
⑥6		DISCONNECT		
⑥7		MOBILE RADIATION SHIELD		
⑥8		MAIN DISCONNECT CONTROL CAT NO. E4502SN		
⑥9		X-RAY ROOM WARNING LIGHT CONTROL PANEL, GE CAT. NO. E4500AK FOR WARNING LIGHT CONTROL ONLY.		

ENVIRONMENT

Ambient Operation Temperature: 68-75° F [20-24° C]
 Allowable Temperature Change: 5° F [3° C] per hour
 Humidity: 40-60 percent non condensing
 Allowable Humidity Change: 5 percent per hour
 Max - Altitude: N/A

ROOM PLAN LOCATED ON PAGE 7-66

ITEM	QUANTITY	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
①	1	XELERIS WORKSTATION	55 lbs	255 btu
②	1	COLOR PRINTER		
③	1	UPS SYSTEM	33 lbs	4436 btu
④	1	VENTRI GANTRY	1432 lbs	2904 btu
⑤	1	VENTRI TABLE	705 lbs	2047 btu
⑥	1	PATIENT POSITIONING MONITOR		
⑦	1	PATIENT POSITIONING CAMERA		
⑧	1	OPERATORS CONSOLE ON CART	99 lbs	273 btu
⑨	1	INTEGRATED POWER SUPPLY	171 lbs	341 btu
⑩	1	COLLIMATOR CABINET	88 lbs	
⑪	1	ACQUISITION COMPUTER		679 btu
***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS *****				
60		TABLE		
61		MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 36 IN. W x 83 IN. H [914mm x 2108mm], CONTINGENT ON A 53 IN. [1346mm] CORRIDOR WIDTH		
62		MAIN DISCONNECT, REFERENCE JUNCTION POINT 'A' ON SHEET E1 FOR DETAILED DESCRIPTION..		
90		OPERATORS CHAIR		

ENVIRONMENT

Ambient Operation Temperature: 55 TO 75 [15 TO 28 C]
 Allowable Temperature Change: 5° F [3° C] per hour
 Humidity: 20-80 percent non condensing
 Allowable Humidity Change: 5 percent per hour
 Max - Altitude: 8000 FT. [2438 M]